

# cardiofax<sup>C</sup>

Electrocardiograph  
ECG-1150

cardiofax continues to evolve



*Fighting Disease with Electronics*

 **NIHON KOHDEN**

# Simple operation and intelligent analysis

## Compact and lightweight: 1.5kg\* and 50 mm thick

The lightweight body guarantees stress-free mobility and the compact size fits anywhere in the clinic or hospital.

\*Weight not including battery and recording paper

## Full 12-lead ECG display

A large 4.8-inch backlit LCD screen shows all 12-lead traces at the same time. In addition to easier waveform viewing, this helps you save recording paper by confirming the recording integrity of all leads before starting recording.

## Highly reliable 12-lead ECG analysis

ECG-1150 includes Nihon Kohden's ECAPS12C interpretation program. ECAPS12C provides simultaneous 12-lead ECG acquisition and analysis with approximately 200 findings and 5 judgements.

## ECG data management


Up to 40 ECG files can be stored in the internal memory.

Approximately 3,000 ECG files can be stored on a 64MB SD card.

With the SD card, you can review and print the stored ECG files.


## AC and battery operation

A rechargeable and ecological NiMH battery provides more than 120 minutes of measurement and recording.

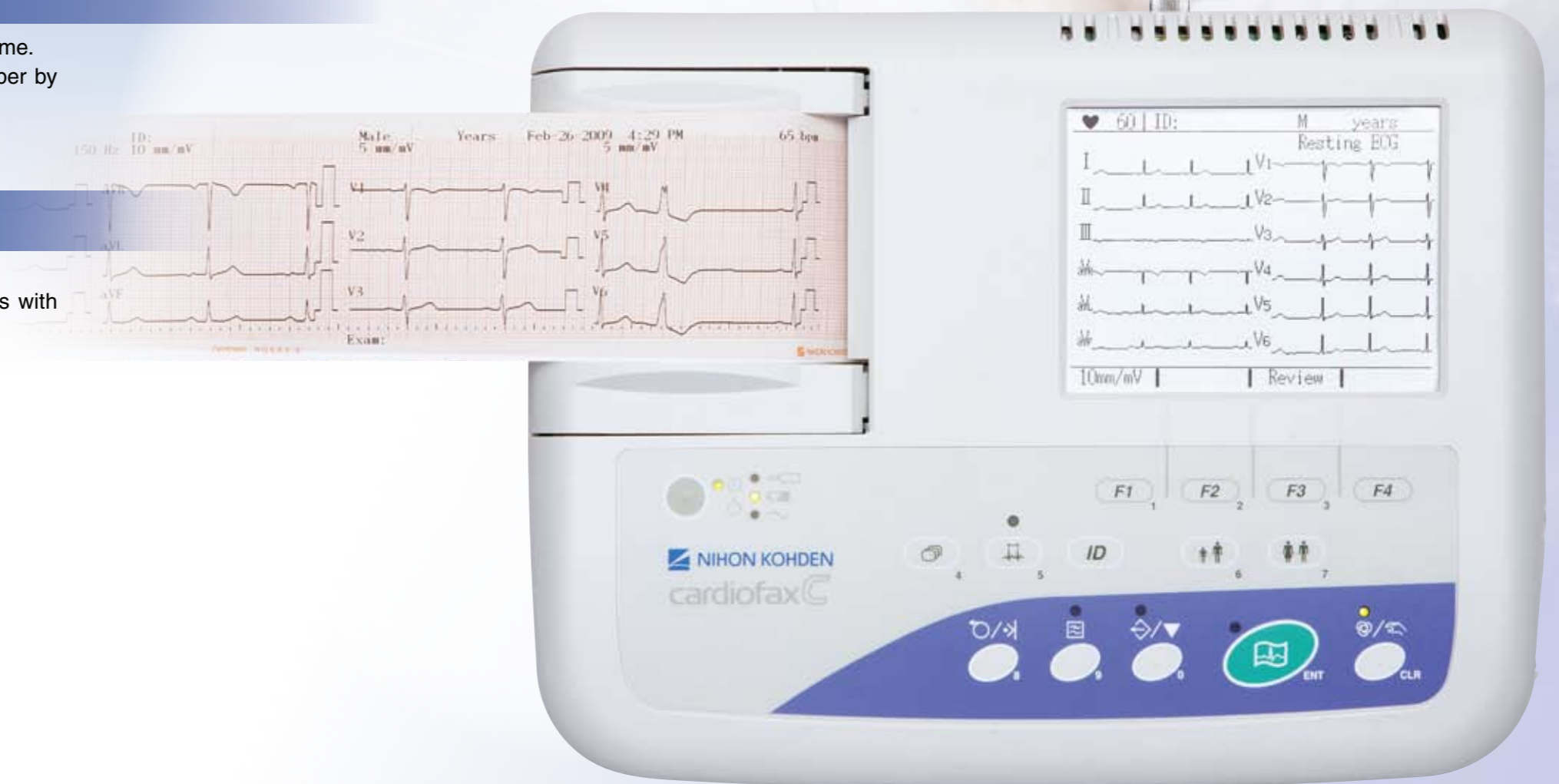



# Recording Examples

## 12 lead ECG (3 traces x 4 groups) with analysis results

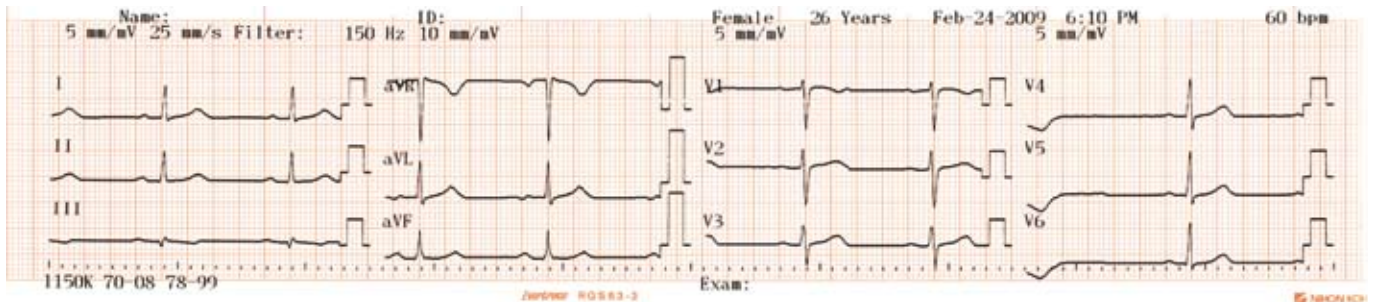

Name:	Feb-24-2009 6:10 PM		1100 Sinus rhythm
ID:	Sex: Female 26 Years		1570 with occasional ventricular premature complexes
Vent rate	60 bpm		9140 ** abnormal rhythm ECG **
PR int	180 ms		
QRS dur	84 ms		
QT/QTc int	400/ 400 ms		
P/QRS/T axis	37/ 30/ 20		
RV5/SV1 amp	1.680/ 1.560 mV		
RV5+SV1 amp	3.240 mV		
Unconfirmed Report			
Reviewed by:			
Exam:			

## 1 channel rhythm recording

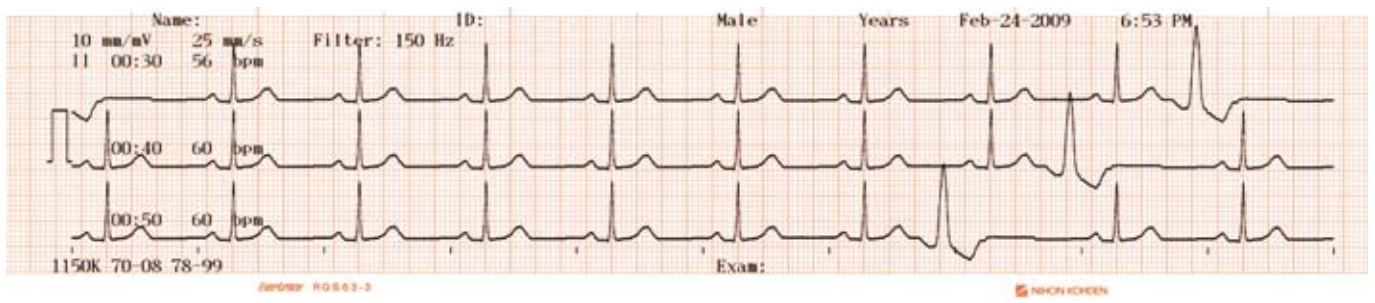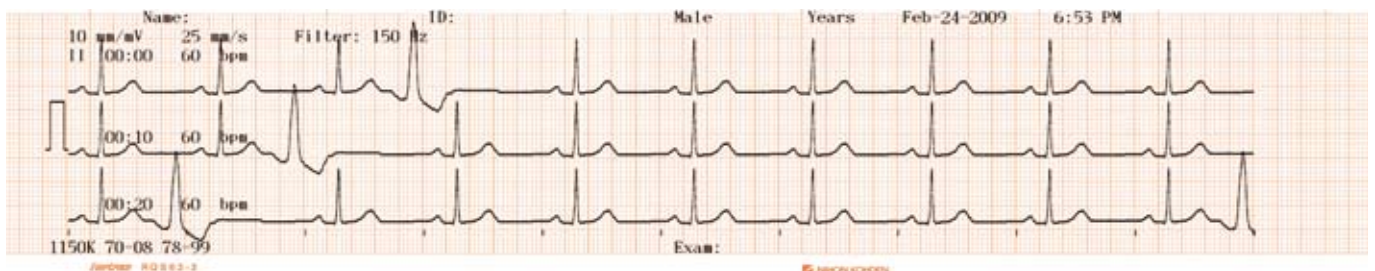

# Options

- Battery
- Cart
- Patient cable hanger
- 63 mm recording paper
- SD card

For more information about accessories and consumables, see Technical Data separately available.



**NIHON KOHDEN CORPORATION**  
 1-31-4 Nishiochiai, Shinjuku-ku, Tokyo 161-8560, Japan  
 Phone +81 (3) 5996-8036 Fax +81 (3) 5996-8100  
[www.nihonkohden.com](http://www.nihonkohden.com)

Printed in Japan on FSC certified paper